
# PRODUCT INFORMATION OTICON DYNAMO SP10, SP8, SP6 AND SP4

## Oticon Dynamo

**Oticon Dynamo** is a BTE 13 super power instrument built on the Inium Sense platform delivering 143dB SPL/82dB full on gain and Inium Sense feedback Shield<sup>5P</sup> an efficient anti-feedback system. Dynamo is designed for severe to profound hearing loss patients who want a slim super power instrument offering easy access to volume control and providing individual listening preferences to ensure optimal speech understanding.

Dynamo features a unique combined approach to signal processesing: Speech Rescue<sup>TM</sup> performing precise frequency composition to ensure full access to high frequency sounds whilst Speech Guard E's advanced adaptive compression preserves the important speech details.

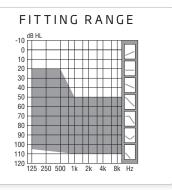

#### Speech Rescue™

Oticon's innovative frequency composition technology, Speech Rescue, increases speech understanding by rescuing speech cues that would otherwise be lost. Speech Rescue positions the high frequency energy in the available space without disturbing the mid frequencies. In this way, it is possible to turn down the inaudible high frequency gain or to give the patient the complete bandwidth of amplification even when Speech Rescue is on.

#### **YouMatic Premium**

YouMatic is a personal automatic system programmed to the client's individual needs and sound preferences. YouMatic controls the sound processing across multiple environments by adjusting the response, directionality, noise management, transient management and compression speed.

#### Speech Guard E

Speech Guard E is the amplitude compression system in in Dynamo SP10 and SP8. It combines two methods of amplification: non-linear and linear in a single compression system that helps preserve the fine details of sound - from soft to loud - making sounds audible, comfortable and clear.

This compression system guards the dynamic contrast of the entire signal to let the rescued high frequency speech sounds appear clear.

#### Inium Sense feedback shield<sup>SP</sup>

Inium Sense feedback shield<sup>sp</sup> ensures that feedback is controlled in the frequency regions where the risk of feedback is highest. It is specifically adapted to severe and profound hearing losses allowing frequency shift across an extended bandwidth.

#### **Family Features**

- Speech Guard E
- O Speech Rescue™
- O Inium Sense feedback shield<sup>sp</sup>
- O Spatial Noise Management
- O Binaural Synchronisation
- O Binaural Coordination
- O YouMatic
- O Dynamic Speech Enhancement (DSE)
- Fitting Bandwidth 6.5 kHz\*
- O Free Focus
- O TriState Noise Management
- O Transient Management
- O Memory (datalogging)
- O Bimodal supportO T-coil
- O 1-COII
- O Special purpose programs
- O AutoPhone Program
- O Back dir
- O Power Bass (streaming)
- O Music Widening (streaming)
- O Multi-band Adaptive Directionality
- O NAL, NAL-RP, DSL
- O 16 Frequency channels
- O ConnectLine, Streamer and remote control
- O DAI input and FM option
- O Voice Priority*i*™
- O In-situ audiometry
- O EasyRECD™
- IP 58 classification: dust and water resistant
- O Analogue volume control w/ mute
- \* Bandwidth accessible for gain adjustments during fitting





People First is our promise to empower people to communicate freely, interact naturally and participate actively









L58474UK / 09.15 Printed on 100% recycled pape

www.oticon.com

#### **PRODUCT OVERVIEW**

Features	Dynamo SP10	Dynamo SP8	Dynamo SP6	Dynamo SP4
Non-linear fitting formulas	DSEsp, NAL, DSL	DSEsp, NAL, DSL	DSEsp, NAL, DSL	DSEsp, NAL, DSL
DSE personal profiles	5	3	3	1
Linear fitting formulas	DSElin, NAL-RP	DSElin, NAL-RP	DSElin, NAL-RP	DSElin, NAL-RP
Speech Rescue™	Yes	Yes	Yes	Yes
Speech Guard E	Yes	Yes	No	No
Single compression	No	No	Yes	Yes
Spatial Noise Management	Yes	No	No	No
Binaural Synchronisation (automatics)	Yes	Yes	No	No
Binaural Coordination (PB operations)	Yes	Yes	Yes	Yes
YouMatic	Premium	Advanced	Essential	No
Noise Management	Transient, Tri state	Tri State	Standard	Standard
Inium Sense feedback shield <sup>sp</sup>	Yes	Yes	Yes	Yes
Free Focus	Premium	Advanced	Essential	Basic
Back dir	Yes	Yes	Yes	No
Power Bass	Yes	No	No	No
Music Widening	Yes	No	No	No
Voice Priority <i>i</i> ™	Yes	Yes	Yes	No
Special Purpose programs (music, lecture etc.)	Yes	Yes	Yes	No
Easy RECD™	Yes	Yes	Yes	No
Fitting Bands	9	8	6	4

#### GENERAL FITTING

Oticon Dynamo instruments are programmed using the Genie 2015.2 fitting software or higher compatible with NOAH 3 or higher.

#### Wireless fitting - FittingLINK

FittingLINK provides a wireless link (Bluetooth) between the PC and one of two wireless enabled hearing instruments, in addition FittingLINK can be used via a USB cable connected to the PC.

#### Cabled Fitting

Use programming cable #3.

#### OPTIONS AND ACCESSORIES

Tamper resistant	Available in 7 colour
battery drawer	

**Damper** Damping element for replacement

**DAI Adaptor** AP 900

FM Receiver

Dedicated

FM Adaptor FM 9

Compatible with Amigo R2 and other universal receivers

Amigo R12 (available in all colours)

Amigo FM T5, T30, T31 Transmitters

#### CONDITIONS

Operating conditions

Temperature: +1°C to +40°C Relative humidity: 5% to 93%, noncondensing

Storage and transportation conditions

Temperature and humidity shall not exceed the below limits for extended periods during transportation and storage: Temperature: -25°C to +60°C

Temperature: -25°C to +60°C Relative humidity: 5% to 93%, non-

condensing

## BTE SUPER POWER OTICON DYNAMO SP10, SP8, SP6 AND SP4



Scale 1:1

**Technical Information**Omnidirectional mode is used unless otherwise stated.

### Warning to the instrument dispenser

Battery size 13

Battery life, calculated\*


Battery life, actual IRIL (IEC 60118-13)


The maximum output capability of the hearing instrument may exceed 132 dB SPL (IEC 711). Special care should be exercised in selecting and fitting the instrument as there may be risk of impairing the remaining hearing of the hearing instrument user.

#### EAR SIMULATOR

Measured according to IEC 60118-0 (1983) and 60711 (1981) and DIN 45605.

OSPL90

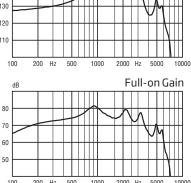

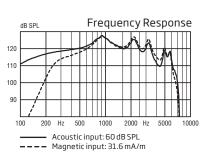



#### 2CC COUPLER


Oticon **Dynamo** 

Measured according to ANSI S3.22 (2003) and S3.7 (1995), IEC 60118-7 (2005) and IEC 60318-5 (2006).

OSPL90

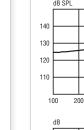

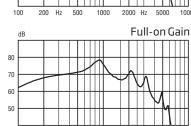

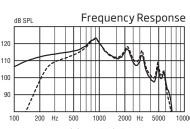

139 dB SPL

129 dB SPL

131 dB SPL

78 dB

67 dB

69 dB

54 dB

100-6200 Hz

107/107 dB SPL 2.0 %

> <2 % 2.0 %

18 dB SPL

33 dB SPL

1.2 mA

1.9 mA

---- Acoustic input: 60 dB SPL
---- Magnetic input: 31.6 mA/m

OSPL90	Peak	143 dB SPL	
	1600 Hz	136 dB SPL	
	Average	138 dB SPL	
Full-on gain	Peak	82 dB	
	1600 Hz	75 dB	
	Average	77 dB	
Reference test gain		60 dB	
Frequency range		100-6500 Hz	
Telecoil output (1600 Hz)	1 mA/m field	104 dB SPL	
1	0 mA/m field	125 dB SPL	
	SPLITS L/R	-	
Total harmonic distortion	500 Hz	6.0 %	
(Input 70 dB SPL)	800 Hz	<2 %	
	1600 Hz	2.0 %	
Equivalent input noise level (A)	) Omni	19 dB SPL	
	Dir	32 dB SPL	
Battery consumption	Quiescent	1.2 mA	
	Typical	1.3 mA	

(IEC PR48)	
* IEC 60118-0 §7.11	246 hours
	80-175 hours
GSM/DECT	800/1400/2000MHz:21/<16/<16 dB SPL

<sup>\*)</sup> Based on the standardised battery consumption measurement (IIC 60118-0.) The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment

